

USER MANUAL

Custom-Built Biaxial Stretcher

Duet 2 WiFi / RepRapFirmware control system

Manual version	1.0
GUI software version	9.6.1
Document date	August 13, 2026
Controller	Duet 2 WiFi / Ethernet family
Primary local interface	Biaxial Stretcher PyQt6 GUI over USB serial
Network interface	Duet Web Control (DWC) in a web browser
GitHub repository	https://github.com/ucerk/stretcher.git

This manual explains how to connect to the stretcher, use the custom graphical user interface, configure Wi-Fi, access Duet Web Control, perform routine movements and heating operations, and shut down the system safely.

Contents

1 Purpose and system overview	2
2 Safety information	2
3 Computer and connection requirements	3
4 Connecting to the stretcher by USB	3
4.1 Physical connection	3
4.2 Selecting the serial port	3
5 Connecting the Duet and computer to the same Wi-Fi network	4
5.1 Option A: phone hotspot	4
5.2 Option B: laboratory or local Wi-Fi network	5
5.3 Re-enabling a previously saved Wi-Fi network	5
5.4 Resetting saved Wi-Fi credentials	5
6 Finding the Duet IP address	5
7 Opening Duet Web Control	6
8 GUI overview	7
8.1 Position and temperature readouts	7
8.2 SD-Macro Mode	7
9 Setup controls	7
10 Heating controls	8
11 Motion controls	8
11.1 Stroke	8
11.2 Continuous separation speed	8
11.3 Timed stepped move	9
11.4 X and Y jog buttons	9
11.5 STRETCH and RETRACT	9
12 Pause, resume, progress, and emergency stop	9
12.1 Pause and resume	9
12.2 Progress bar	9
12.3 Hard emergency stop	10
13 Calibration control	10
14 Console	10
15 Recommended normal operating sequence	11
16 Using Duet Web Control	11
16.1 The main DWC pages	12
16.2 User macros	12
16.3 Creating a new macro in DWC	13
16.4 Editing an existing macro	13
16.5 Running a macro	14
16.6 Changing config.g	15

16.7 Using the DWC Console 15

17 Troubleshooting 16

18 Useful Duet commands 17

19 Shutdown and storage 17

20 Software and configuration notes 17

21 Reference documentation 18

1. Purpose and system overview

The biaxial stretcher is controlled by a Duet controller running RepRapFirmware. The custom desktop GUI communicates with the Duet through a USB serial connection and provides simplified controls for motion, heating, calibration, Wi-Fi setup, pause/resume, and emergency stopping.

The system can be accessed in three complementary ways:

1. **Custom GUI over USB** – the normal interface for day-to-day operation of the stretcher.
2. **Duet Web Control (DWC)** – a browser-based interface available when the Duet and the computer are connected to the same network.
3. **PanelDue** – if installed, provides local status and selected job-control functions such as pause/cancel for SD-card jobs.

Note: The USB GUI and Duet Web Control may both be connected at the same time. However, do not issue competing movement commands from two interfaces simultaneously. Use one interface as the active control point for motion.

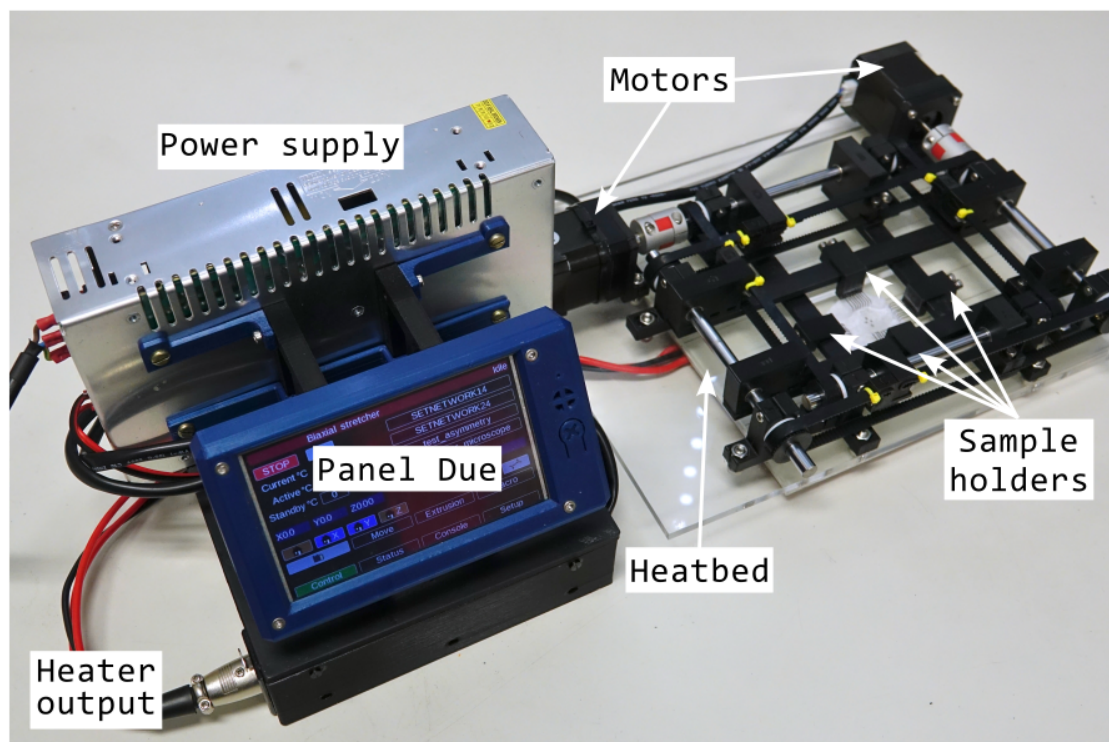


Figure 1: Complete biaxial stretcher, Duet controller (hidden), and Panel Due.

2. Safety information

WARNING: The stretcher contains powered motors, moving grips, stored elastic energy, a heater, and exposed electronics. Keep hands, tools, cables, and loose material clear of moving parts whenever the motors are enabled.

Before each use:

- inspect the frame, grips, leadscrews, couplings, cables, heater wiring, and controller;
- confirm that no object can obstruct X or Y motion;

- confirm that the specimen and fixtures are mounted securely;
- verify that the selected stroke and speed are appropriate before pressing a movement button;
- verify that the temperature setpoint is suitable for the experiment;
- know how to remove main power if the software emergency stop is insufficient.

The GUI contains a **HARD EMERGENCY STOP (M112)** button. This sends RepRap-Firmware command M112. After an emergency stop, the Duet must be reset before normal operation can continue. The GUI indicates that the PanelDue STOP/reset procedure must be completed before reconnecting.

3. Computer and connection requirements

For normal operation you need:

- the biaxial stretcher and Duet controller;
- a USB cable between the Duet and the computer;
- the Biaxial Stretcher GUI;
- for network access, a 2.4 GHz Wi-Fi network or phone hotspot;
- a modern web browser for Duet Web Control.

The Duet 2 WiFi uses 2.4 GHz Wi-Fi. A 5 GHz-only hotspot or wireless network cannot be used by the Duet. If a phone hotspot is used, enable a 2.4 GHz or compatibility mode if the phone provides that option.

4. Connecting to the stretcher by USB

4.1. Physical connection

1. Make sure the stretcher is in a safe mechanical state and that no movement will cause a collision.
2. Connect the Duet to the computer using the USB cable.
3. Apply the normal stretcher power supply when motor or heater operation is required.
4. Start the Biaxial Stretcher GUI.

4.2. Selecting the serial port

The GUI automatically lists serial ports detected by `pyserial`.

Note: If the GUI is started before connecting the stretcher to the computer/laptop, the serial ports won't be listed! Close and restart the software.

Typical names are:

- **Windows:** COM3, COM4, etc.
- **macOS:** `/dev/tty.usbmodem...`


Select the Duet port from the drop-down list and press **CONNECT**.

When the GUI connects, it sends initial controller commands that:

- place Wi-Fi in idle mode with M552 S0;
- allow movement with M564 H0 S0;
- enable the motors with M17.

WARNING: Because pressing **CONNECT** enables the motors, connect only when the mechanism is clear and safe to energize. Also note that the motors may jolt upon energizing which could affect the sample!

Note: The GUI intentionally sends M552 S0 when connecting. Therefore, a Wi-Fi connection that was previously active may be placed into idle mode. Re-enable Wi-Fi using the GUI after the USB connection is established.

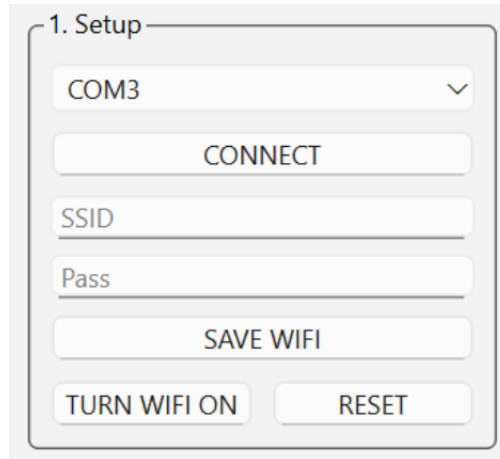

The screenshot shows a GUI window titled "1. Setup". It contains a dropdown menu set to "COM3" with a checkmark icon. Below this is a "CONNECT" button. Further down are two input fields labeled "SSID" and "Pass". Below these fields is a "SAVE WIFI" button. At the bottom of the window are two buttons: "TURN WIFI ON" and "RESET".

Figure 2: Setup module of the GUI with controls for connecting to the Duet Due (via a USB cable) and WiFi (for configuring the Duet).

5. Connecting the Duet and computer to the same Wi-Fi network

5.1. Option A: phone hotspot

A phone hotspot is convenient when no laboratory router is available or when the laboratory network cannot be used. Institutional networks may require device registration, browser-based authentication, user credentials, or other security settings that are not supported by the Duet. In addition, network segmentation or wireless-client isolation may prevent the laptop from accessing the Duet even when both devices appear to be connected. In these situations, a private phone hotspot provides a simpler local network for the laptop and Duet.

Note: An unfamiliar IP address is not necessarily an error. Laboratory networks may use a different private address range from the expected 192.168.x.x range. However, if the Duet does not obtain an address, or if the laptop and Duet are placed on isolated network segments, Duet Web Control may not be accessible. Use a private phone hotspot instead.

1. Enable the phone hotspot.
2. Ensure that the hotspot is operating in a mode compatible with 2.4 GHz Wi-Fi.
3. Connect the laptop to the phone hotspot in the normal way.
4. Connect the laptop to the Duet by USB and open the Biaxial Stretcher GUI.
5. Press **CONNECT** in the GUI.
6. Enter the hotspot name into the **SSID** field.
7. Enter the hotspot password into the **Pass** field.

8. Press **SAVE WIFI**.
9. Allow several seconds for the Duet to join the hotspot.

The GUI's **SAVE WIFI** function performs the equivalent of:

```
M552 S0
M587 S"network-name" P"network-password"
M552 S1
```

The M587 command stores the network credentials in the Duet Wi-Fi module. The credentials therefore do not normally need to be entered after every power cycle.

5.2. Option B: laboratory or local Wi-Fi network

The procedure is the same as for a phone hotspot:

1. Connect the laptop to the required local Wi-Fi network.
2. Confirm that the network provides 2.4 GHz Wi-Fi access.
3. Connect to the Duet through USB in the GUI.
4. Enter the local network SSID and password.
5. Press **SAVE WIFI**.

5.3. Re-enabling a previously saved Wi-Fi network

If the credentials are already stored, it is not necessary to press **SAVE WIFI** again. After connecting by USB, press the GUI button labelled **WIFI: OFF**. The button sends M552 S1 and changes to **WIFI: ON**.

5.4. Resetting saved Wi-Fi credentials

The GUI **RESET** button beside the Wi-Fi control clears all stored Wi-Fi network credentials. It first places the Wi-Fi module in idle mode, allowing the reset command to be processed, and then disables the module after the reset. It sends commands equivalent to:

```
M552 S0
G4 P500
M588 S"*"
G4 P500
M552 S-1
```

WARNING: The Wi-Fi RESET button erases remembered Wi-Fi networks. Use it only when you intentionally want to remove the stored SSID/password configuration.

6. Finding the Duet IP address

Duet Web Control is opened using the Duet's network address. When DHCP is used, the router or phone hotspot assigns that address automatically.

Possible ways to find the IP address are:

- note the IP address reported by the Duet when Wi-Fi connects. It appears in the GUI console or in the PanelDue Console;
- view the connected-device list in the router or phone hotspot;
- connect by USB with a serial terminal and send M552.

Sending M552 with no parameters reports the current network state and IP address.

Note: The custom GUI does not currently provide a general-purpose G-code command entry box. If you need to send M552 manually over USB, close the GUI first so another serial terminal, for example YAT, can open the Duet serial port.

7. Opening Duet Web Control

Once the laptop and Duet are on the same network:

1. Open a web browser on the laptop.
2. Enter the Duet IP address, for example `http://192.168.1.90/`.
3. Alternatively, try `http://Stretcher.local/`.
4. Duet Web Control should load in the browser.

Duet Web Control provides access to controller status, files, configuration, the G-code console, heater controls, movement controls, job controls, and firmware/configuration tools.

WARNING: On a standalone Duet, Duet Web Control normally uses HTTP rather than HTTPS. Use it only on a trusted local network or private hotspot. Do not expose the Duet directly to an untrusted or public network.

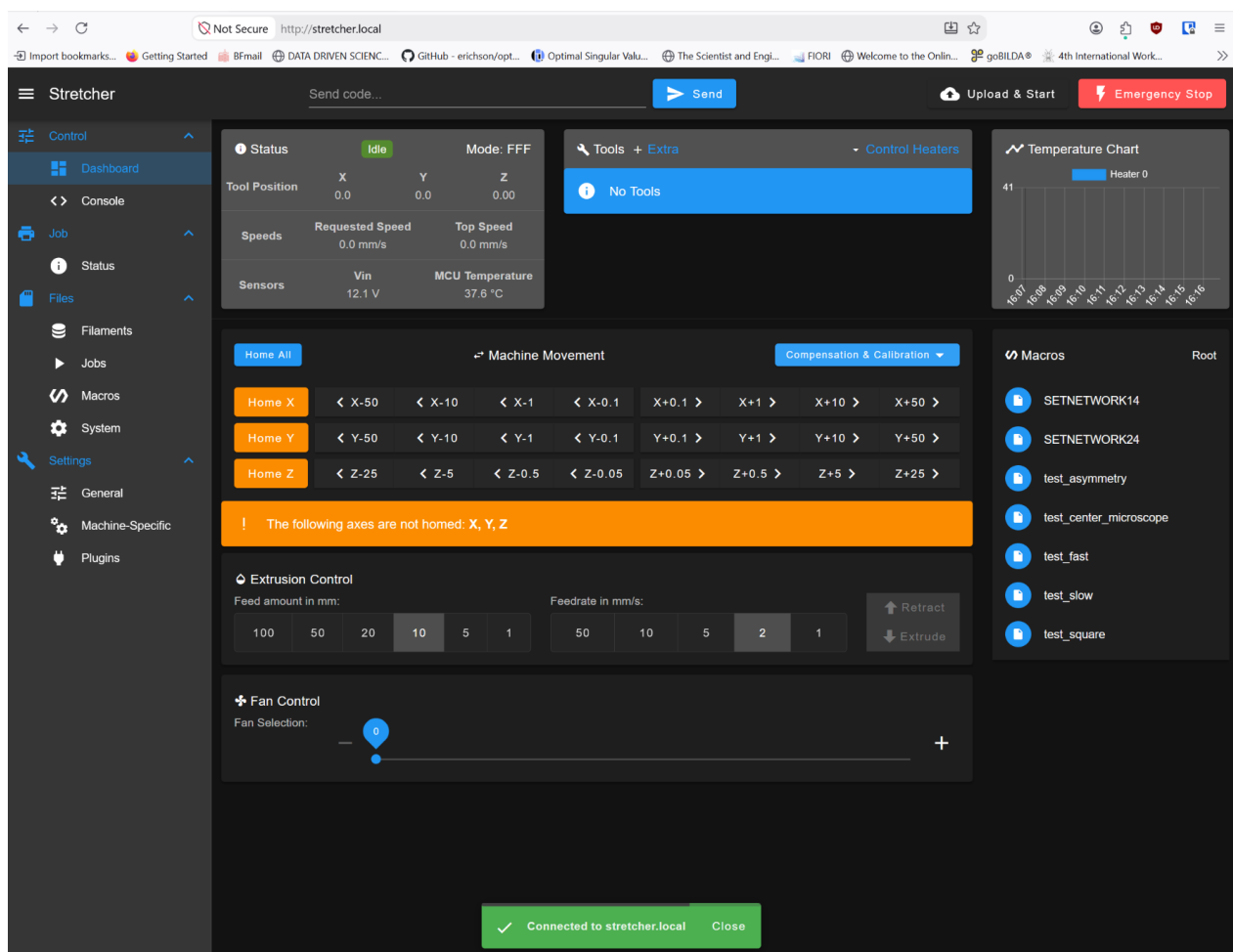


Figure 3: Web based Duet control center.

8. GUI overview

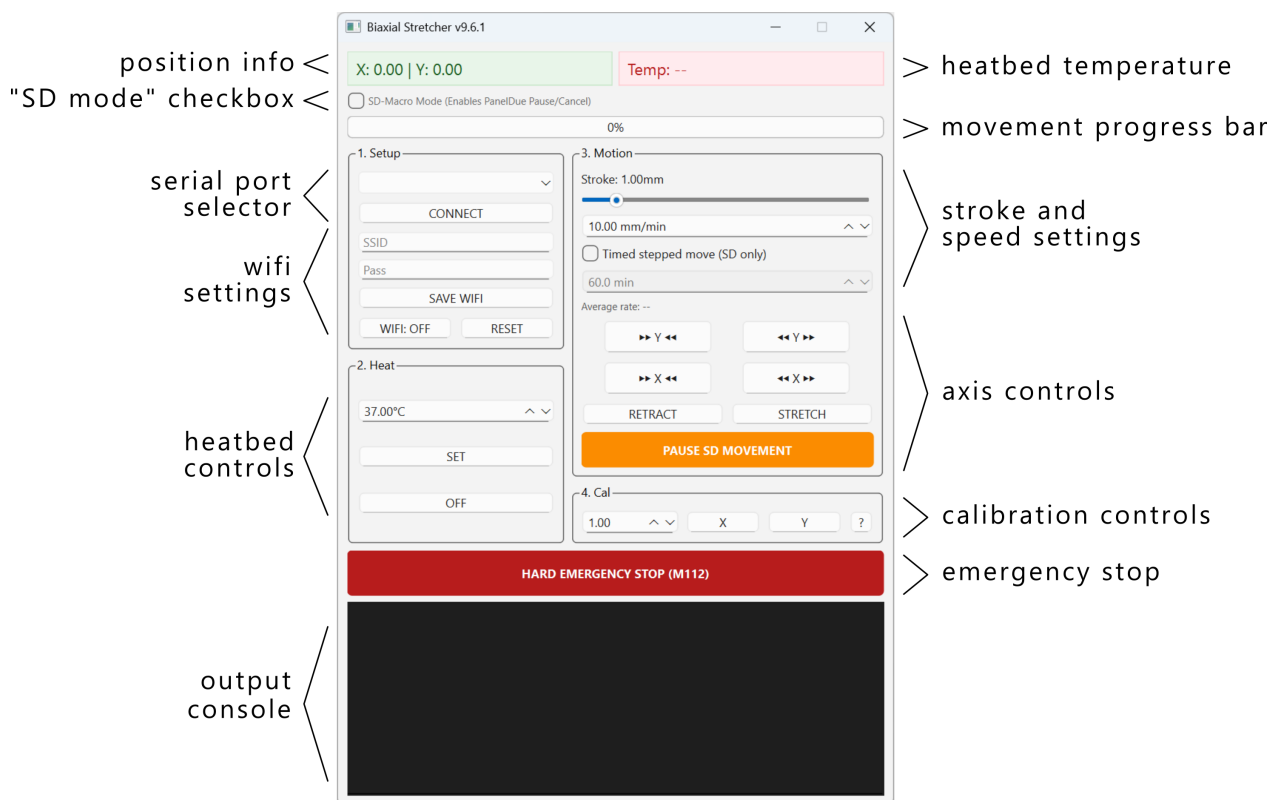


Figure 4: GUI overview.

The GUI is divided into several functional areas.

8.1. Position and temperature readouts

At the top of the window:

- **X / Y** displays the latest coordinates reported by the Duet;
- **Temp** displays the bed/heater temperature reported by the controller.

The GUI periodically requests position, temperature, and controller job status from the Duet.

8.2. SD-Macro Mode

The checkbox **SD-Macro Mode (Enables PanelDue Pause/Cancel)** changes how normal movements are executed.

When enabled, the GUI writes the movement to an SD-card G-code file and starts it as a Duet job. This enables cleaner pause/resume interaction using M25/M24 and PanelDue job controls.

For a simple direct USB move, leave SD-Macro Mode unchecked unless SD pause/cancel behavior is required.

9. Setup controls

The **1. Setup** group contains:

- | | |
|-------------------------|--|
| Serial-port list | Selects the USB serial connection to the Duet. |
| CONNECT | Opens the serial connection and performs the initial controller setup. |

SSID	Wi-Fi network name.
Pass	Wi-Fi password; displayed as hidden characters.
SAVE WIFI	Stores the selected network on the Duet and enables Wi-Fi.
WIFI: OFF / ON	Enables or idles Wi-Fi without changing the saved credentials.
RESET	Erases the Duet's remembered Wi-Fi networks.

10. Heating controls

The **2. Heat** group contains a temperature input and two buttons.

- The temperature control range is 0 °C to 40 °C.
- The default displayed setpoint is 37 °C.
- Pressing **SET** sends M140 H0 S<temperature>.
- Pressing **OFF** sends M140 H-1 to disable the heater.

The temperature displayed at the top of the GUI is the temperature reported by the Duet. The Python program does not perform PID control itself; temperature regulation is handled by RepRapFirmware on the Duet.

WARNING: Confirm the correct heater and sensor configuration before use. Do not rely on the GUI setpoint alone when validating a temperature-sensitive experiment.

11. Motion controls

11.1. Stroke

The horizontal **Stroke** slider defines the commanded change in separation.

The current software uses a slider range of 1 to 1000 internal units, with 100 units corresponding to 1.00 mm. Therefore the selectable displayed range is approximately 0.01 mm to 10.00 mm.

For the slider keyboard controls:

- ←/→: decrease/increase by 0.01 mm;
- ↓/↑: decrease/increase by 0.10 mm.

11.2. Continuous separation speed

The speed field defines the requested **separation speed** in mm/min. The software converts this to the required motor feed rate internally.

The current continuous-mode minimum is 1.20 mm/min. Lower average rates must use **Timed stepped move** mode.

Note: In the currently supplied v9.6.1 source, the custom speed-box keyboard handler uses Up/Down increments of 1.0 and Left/Right increments of 0.5. If the intended behavior is Up/Down = 0.10 and Left/Right = 0.01, update the source code available on github.

11.3. Timed stepped move

Enable **Timed stepped move (SD only)** when an average separation rate below the continuous RepRapFirmware motion floor is required.

WARNING: The time stepped motion is **not continuous**, but happens in internally calculated steps based on the calibration.

When enabled:

- SD mode is forced on;
- the normal speed box is disabled;
- the duration field becomes active;
- the GUI distributes complete motor microsteps across the selected duration;
- the average rate is displayed as stroke divided by duration.

The duration control supports approximately 1 to 1440 minutes.

11.4. X and Y jog buttons

The four jog buttons move either the X or Y axis using the currently selected stroke and movement mode. The arrow symbols indicate whether the grips move in the positive or negative configured direction.

Perform the first jog after setup at a small stroke and verify the physical direction before using a larger movement.

11.5. STRETCH and RETRACT

STRETCH and **RETRACT** command both X and Y axes together.

For a biaxial move, the software compensates the vector feed rate so that the individual X and Y motor speeds correspond to the requested separation speed.

WARNING: Always verify the selected stroke and direction before pressing **STRETCH** or **RETRACT**. Movement begins when the command is sent; there is no additional confirmation dialog.

12. Pause, resume, progress, and emergency stop

12.1. Pause and resume

The **PAUSE SD MOVEMENT** button is enabled only during SD-card movements.

- Pause sends M25.
- Resume sends M24.
- The GUI monitors the Duet job state and updates the button as the controller enters pausing, paused, resuming, or processing states.

Direct USB movements cannot use the same SD-job pause mechanism.

12.2. Progress bar

For normal moves, the GUI estimates progress from the requested distance and speed. For timed stepped protocols, progress is reported by messages generated by the SD job itself.

12.3. Hard emergency stop

Pressing **HARD EMERGENCY STOP (M112)** immediately sends the Duet emergency-stop command and closes the serial connection.

After an emergency stop:

1. ensure that all physical motion has stopped;
2. remove the cause of the emergency;
3. reset the Duet using by pressing the **STOP** button in the top-left corner of the Panel Due;
4. reconnect only when the mechanism is safe.

WARNING: After an emergency stop, reconnecting to the device re-energizes the stepper motors. The motors may make a **small, sudden movement** as the rotor aligns with the magnetic field established by the stepper driver. This movement may affect the sample.

13. Calibration control

The **4. Cal** group provides a simple displacement calibration adjustment.

1. Command a known stroke.
2. Measure the actual physical separation change using a suitable measurement tool.
3. Enter the measured value in the calibration input box.
4. Press **X** or **Y** for the axis being calibrated.

The software calculates a scale factor:

$$\text{scale} = \frac{\text{commanded stroke}}{\text{measured stroke}}$$

and sends an **M92** command that scales the **current live** steps-per-mm value for the selected axis.

Note: The calibration takes effect immediately but is not automatically saved to the permanent Duet configuration. To retain the calibrated steps-per-mm value after a restart, update the corresponding **M92** value in **config.g** via Duet Web Control; see [Opening Duet Web Control](#).

14. Console

The black console area at the bottom of the GUI shows status messages, network information, calibration reports, movement events, and errors. Repetitive position, temperature, and **ok** responses are filtered to keep the display readable.

Use the console when diagnosing:

- serial connection failures;
- Wi-Fi connection messages;
- SD upload errors;
- timed-protocol errors;
- calibration responses;
- controller warnings.

15. Recommended normal operating sequence

1. Inspect the stretcher and clear the motion area.
2. Connect the Duet to the computer by USB.
3. Connect to power supply (to operate the motors/heater).
4. Start the GUI.
5. Select the Duet serial port and press **CONNECT**.
6. If network access is required, enable Wi-Fi or save the current network credentials.
7. Verify the live position and temperature readouts.
8. If required, set the heater temperature and allow the system to stabilize.
9. Select a small stroke and verify X/Y motion direction before mounting or loading a critical specimen.
10. Mount the specimen using the approved fixture procedure.
11. Select stroke, speed, SD mode, or timed stepped mode as required.
12. Execute the required X/Y or biaxial movement.
13. Monitor the specimen, temperature, position, progress bar, and controller messages throughout the movement.
14. Use SD pause/resume (only for SD-card movements).
15. At the end of the session, switch the heater off and return the mechanism to a safe unloaded state.
16. Close the GUI normally. On close, the program requests M18, disabling the X/Y motors and removing holding current before the serial worker exits.
17. Disconnect the stretcher power supply.

16. Using Duet Web Control

Duet Web Control (DWC) provides direct access to the Duet controller through a web browser. It is useful for controller-level diagnostics, configuration, file management, and operations that are intentionally not exposed in the simplified stretcher GUI.

Typical uses include:

- viewing controller state, positions, and temperatures;
- sending diagnostic G-code commands from the **Console**;
- checking network state using M552;
- inspecting or editing `config.g`;
- viewing and managing files stored on the Duet SD card;
- creating, editing, and running macros;
- reviewing heater, motor, and machine configuration;
- checking firmware and controller information.

WARNING: DWC provides lower-level control than the custom GUI. Commands entered in the Console and changes made to configuration or macro files can cause immediate motor movement, heater operation, or changes to controller behaviour. Only trained users should make such changes.

16.1. The main DWC pages

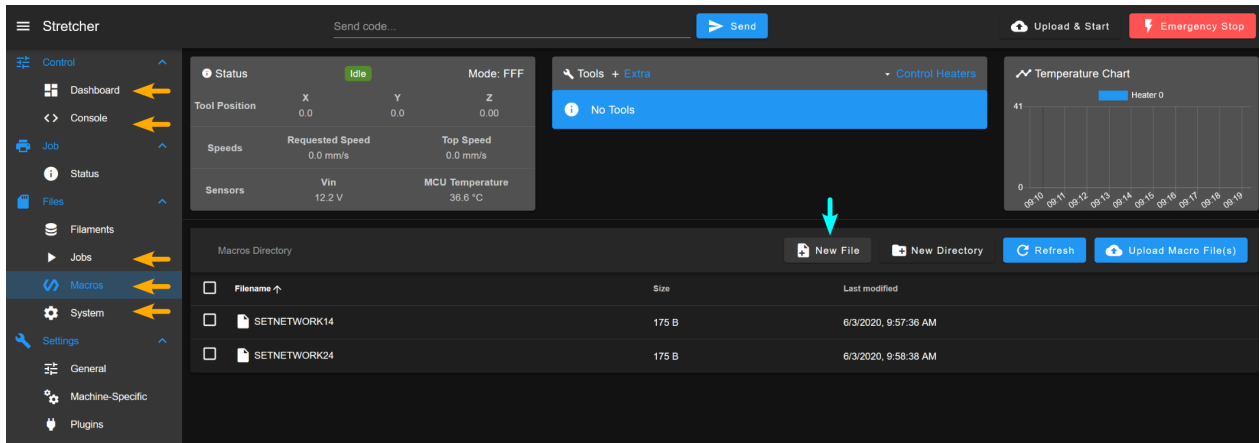


Figure 5: DWC main pages.

The exact appearance of Duet Web Control depends on the DWC version, but the most useful areas for operating and maintaining the stretcher are:

- Dashboard** Shows the current machine state and provides controller-level movement and heater controls.
- Console** Allows G-code commands to be sent directly to RepRapFirmware and displays the controller response.
- Macros** Displays user-defined macro files stored in the Duet's `/macros/` directory. Macros can be created, edited, uploaded, organised, and executed from this page.
- Jobs / Files** Provides access to G-code jobs and other files stored on the SD card.
- System** Provides access to system configuration files stored in the `/sys/` directory, including `config.g`.

16.2. User macros

A macro is a text file containing one or more G-code commands that can be executed as a sequence. Macros are useful when the same controller operation needs to be performed repeatedly.

Normal user macros are stored in:

`/macros/`

The filename is used as the macro name displayed in Duet Web Control and **on PanelDue**. Subdirectories may be used to organise a larger number of macros.

Note: Do not confuse normal user macros in `/macros/` with system files in `/sys/`. Files in `/sys/`, such as `config.g`, homing files, pause/resume files, and other system macros, can affect fundamental machine behaviour.

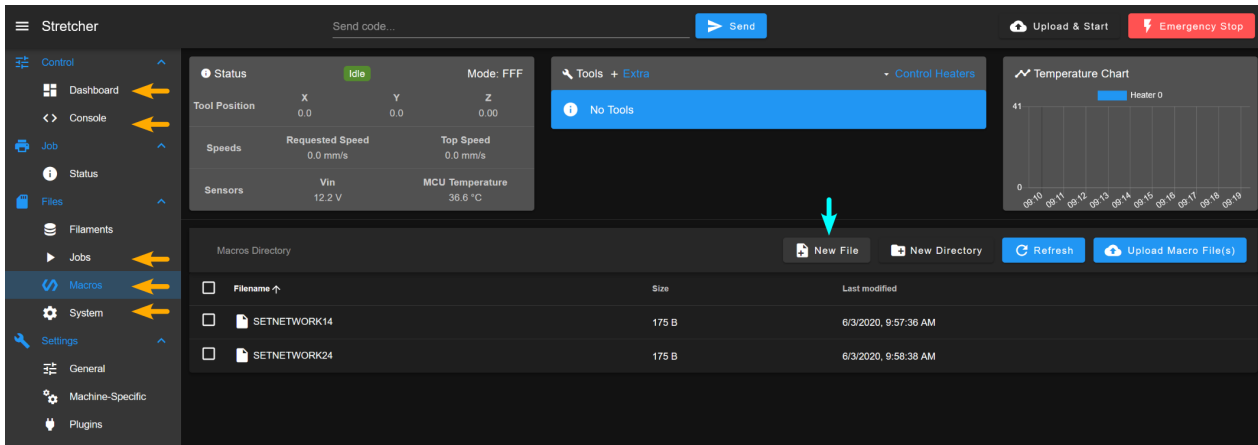


Figure 6: DWC main pages/areas (orange arrows) and how to create a new (macro) file (cyan arrow).

16.3. Creating a new macro in DWC

To create a new macro:

1. Open Duet Web Control as described in [Opening Duet Web Control](#).
2. Open the **Macros** page.
3. Select **New File** (Fig. ??).
4. Enter a descriptive filename for the macro.
5. Enter the required G-code commands into the editor.
6. Save the file.
7. Confirm that the new macro appears in the macro list.

For example, a diagnostic macro named **Stretcher Status.g** could contain:

```
; Stretcher Status
; Report current position
M114

; Report heater temperature
M105

; Report network state
M552
```

Lines beginning with a semicolon (;) are comments and can be used to document what the macro is intended to do.

Note: When creating a new macro, begin with commands that do not cause motion or heating whenever possible. Verify the macro from the Console output before adding commands that actuate the machine.

16.4. Editing an existing macro

Existing user macros can be edited directly from Duet Web Control.

1. Open the **Macros** page.

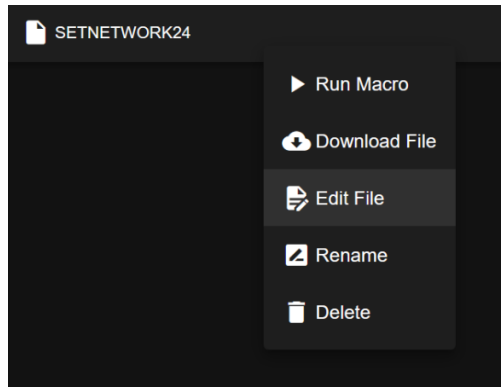


Figure 7: Right click on existing file opens up the drop down menu.

2. Locate the required macro.
3. Open the macro's context menu using R-mouse-click and select **Edit File** (Fig. 7).
4. Review the existing commands before making changes.
5. Modify the required G-code.
6. Save the file.
7. Re-open the file if necessary and verify that the intended changes were saved correctly.
8. Test the macro with the stretcher in a safe mechanical state.

WARNING: Saving a macro changes the copy stored on the Duet SD card. There is no automatic check that a movement command is mechanically safe for the stretcher. Carefully check axis letters, movement direction, distance, feed rate, and heater commands before executing an edited macro.

For important macros, it is recommended to keep a known-good backup before making changes. A macro can be downloaded from DWC and retained on the computer before the SD-card version is edited.

16.5. Running a macro

To execute a user macro, there are two normal options:

1. Run the macro from Duet Web Control

- (a) Open the **Macros** page in DWC.
- (b) Locate the required macro.
- (c) Select the macro.
- (d) R-click and select **Run Macro** (Fig. 7).

2. Run the macro from PanelDue

- (a) Open the **Macros** menu on the PanelDue.
- (b) Locate the required macro.
- (c) Select the macro by pressing it on the screen.

Note: The same macro file is used whether it is started from Duet Web Control or from PanelDue. Changes made to a macro in DWC therefore also affect the version available from PanelDue.

A macro can also be called from another G-code file or from the DWC Console using **M98**. For example:

```
M98 P"/macros/Stretcher Status.g"
```

This executes the commands stored in the specified macro file.

WARNING: Do not run an unfamiliar macro simply to determine what it does. Open the file and inspect its contents first. A macro may contain movement, heater, configuration, or emergency-control commands.

16.6. Changing `config.g`

The main RepRapFirmware machine configuration is stored in `/sys/config.g`. It can be viewed and edited from the **System** page in Duet Web Control.

One relevant example for the stretcher is calibration. The GUI calibration function changes the live X or Y steps-per-mm value using **M92**, but the change is not automatically written to the permanent configuration. After confirming the calibrated value, the corresponding **M92** command in `config.g` can be updated so that the value is restored when the Duet starts.

Before changing `config.g`:

1. Open the **System** page.
2. Download or otherwise record the current `config.g` as a backup.
3. Open `config.g` in the DWC editor.
4. Locate the command that requires modification.
5. Change only the required parameter.
6. Save the file.
7. Restart the Duet when appropriate and verify the resulting configuration before operating the stretcher normally.

For example, if calibration establishes that the required X-axis value is **X1600**, only the appropriate parameter in the existing **M92** line should be changed. Do not copy this example value unless it is the value established for the actual machine.

WARNING: An incorrect `config.g` can change motor scaling, directions, current limits, heaters, sensors, networking, and other fundamental controller settings. Always retain a known-good copy before editing the configuration.

16.7. Using the DWC Console

The DWC **Console** is useful for testing individual G-code commands before placing them into a macro.

For example:

```
M552
```

reports the current network state, while:

M92

reports the currently configured steps-per-mm values.

Testing individual diagnostic commands in the Console can help distinguish a problem with a particular G-code command from a problem with the macro file itself.

WARNING: Commands entered in the DWC Console are sent directly to the controller. Movement or heater commands may take effect immediately after they are submitted.

17. Troubleshooting

Problem	Likely cause	Action
No serial port appears	USB cable, driver, port, or power issue	Reconnect USB, try another data-capable cable/port, and check the operating-system device list.
GUI reports OFFLINE	Wrong serial port or port already in use	Select the correct port. Close other serial-terminal software and reconnect.
Wi-Fi will not connect	Wrong SSID/password, 5 GHz-only network, weak signal	Confirm exact credentials, use 2.4 GHz Wi-Fi, simplify hotspot settings, and retry SAVE WIFI.
Wi-Fi worked before CONNECT but DWC disappears	GUI CONNECT sends M552 S0	Press WIFI: OFF/ON to re-enable Wi-Fi, or use SAVE WIFI if credentials must be changed.
Cannot find DWC address	DHCP assigned an unknown IP	Check router/hotspot clients, use the mDNS name, or query M552 through a serial terminal.
Browser changes HTTP to HTTPS	Browser security behavior	Enter the explicit <code>http://</code> address. Standalone DWC does not normally provide HTTPS.
Axis moves in unexpected direction	Direction/configuration mismatch	Stop immediately and verify axis configuration before continuing.
Timed protocol reports an error	Stroke below a motor step or selected duration too short	Increase the stroke or duration and retry.
Pause button is disabled	Movement is direct USB, not an SD job	Use SD-Macro Mode if pause/resume capability is required.
Temperature differs from target	Heater model, PID tuning, sensor placement, or calibration issue	Verify heater configuration and sensor accuracy; perform controller-level heater tuning where appropriate.
Emergency stop was triggered	Duet is halted	Make the mechanism safe, reset the Duet, then reconnect.

18. Useful Duet commands

Command	Purpose
M552	Report current network state and IP address.
M552 S0	Start/leave Wi-Fi module in idle mode.
M552 S1	Enable Wi-Fi client mode and connect to a remembered network.
M587 S"SSID" P"password"	Store Wi-Fi network credentials.
M587	List remembered Wi-Fi networks.
M588 S"*"	Forget remembered Wi-Fi networks.
M115	Report firmware information.
M140 H0 S37	Example: set heater 0 target to 37°C.
M140 H-1	Disable the configured bed/heater output used by the GUI.
M25	Pause an SD job.
M24	Resume an SD job.
M112	Emergency stop.
M18	Disable motors/remove holding current.

19. Shutdown and storage

1. Finish or abort any active movement.
2. Return the specimen and mechanism to a safe unloaded state.
3. Switch the heater off.
4. Save any required experimental records.
5. Close the GUI normally so the application can send M18 and shut down the serial worker.
6. Switch off the main stretcher power supply.
7. Disconnect USB if the controller is not required to remain powered.
8. Clean the specimen area and inspect for fluid ingress, loose fixtures, or damage.

20. Software and configuration notes

This manual describes GUI software version 9.6.1. If the GUI is changed, review this manual whenever the change affects:

- button functions or labels;
- motion direction, stroke, or speed behavior;
- timed stepped movement;
- heater commands;
- Wi-Fi connection behavior;
- emergency-stop or shutdown behavior;
- calibration calculations.

21. Reference documentation

Official Duet3D documentation used for controller/network procedures:

- [Getting connected to your Duet](#)
- [Setting up networking on Duet](#)
- [Duet Web Control Manual](#)
- [RepRapFirmware G-code dictionary](#)